
Product Specifications

Name: MEJES VII E-ion Molecule removes Speckle Cosmetic Device
Model: MEJES2023FR012388
Charging Type: USB Charging
Power Voltage: 5V~2A
Power Output: 8W
Battery Capacity: 2000mAh
Charging Duration: 3 hours
Battery Working Time: 90 Minutes
Standby Time: 90 Days (If the unit hasn't been used for more than 90 days, please make it fully charged before use)
Size: 218*96*67mm

Product Structures

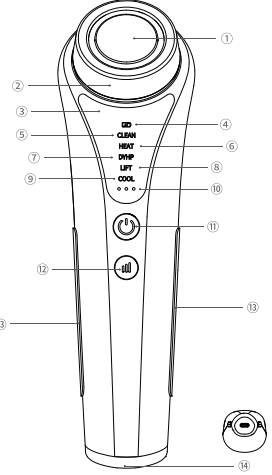① Massager

② Circclip

③ LED

④ Battery

⑤ CLEAN

⑥ HEAT

⑦ DYHP

⑧ LIFT

⑨ COOL

⑩ Gear

⑪ Power/Mode

⑫ Gear

⑬ Sensing

⑭ TYPE-C

Use

Alert

Do not use the device continuously on the same area for a long time, which may cause potential skin problems.

Do not contact lenses around eye areas, which may lead to eye skin problems.

1.Clean the face and keep face moist.

2.Long press power to start the device. Switch to correct mode.

3.Take the device on your face and move gently to massage skin.


Use of neck freckles (Max: 5 min)


Use of facial freckles (Max: 30 min)

1.Along the direction of the arrow, move gently from the chin.

2.Along the chin to the angulus oris and the nasal side to the external canthus move.

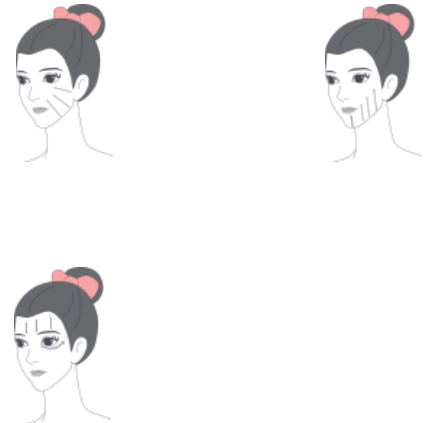Use of eye freckles (Max: 2 min)

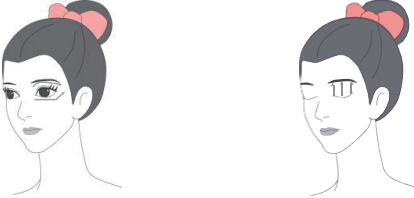FAQ

Press and hold for 2 seconds to turn on/off.

Press the L key to adjust the gear.

Note: It is recommended to use the lowest level for sensitive skin.

MODE

LIFT(E-ion Molecule removes Speckle):Moisten the skin with negative E-Water(Standard Configuration),Enable LIFT (red) mode (Ⓓ),The massage panel gently touches the skin,Slowly moving on the skin;Let positive and negative charge ions collide,Thereby disintegrating melanin.

Important:

1. Keep skin moist with conductive gel when LIFT Mode.

2. After LIFT Mode, cleaning is not allowed within 4 hours

3.Negative E-Water Usage order: 20c~100c~500c~1000c... gradually

Safety Warnings

●Please read all instructions carefully for safety and proper use.

●Use only as directed.

●Immediately stop using if any malfunction or unusual reaction to treatment occurs.

●Precautions listed in this manual are very important. For your safety, please follow them carefully.

●Precautions are classified as WARNING or CAUTION in order to clearly specify the level of damage and possible injury or damage that may occur if the device is misused.

WARNING & CAUTION

Prohibited

Improper use of this device may result in death or injury

Remove Power Cord

Hold the head of the power cord instead of any part of the cord when unplugging the power cord from USB interface. Damage to the cord may cause electrical shock, short circuit or fire hazard.

Always unplug the power cord from USB interface when not use. Do not handle the AC adapter with wet hands. Which may result in electrical shock, short circuit or fire hazard.

If any malfunction occurs, unplug the power cord from the device immediately. Possible malfunctions may include button failure, smoke, unusual noise or smell, please stop using it.

Prohibited

Do not use if the power cord is damaged or if the wall socket is loose which may cause electrical shock, short circuit.

Disabled for children under 17 years of age,Pregnant women are disabled.If you are sick, please follow the doctor's instructions.

Do not use extension cords. This may cause electrical shock, overheating or fire hazard.

Prohibited

Not for use by children, elderly, or disabled person. Improper or unattended use may cause injury.

Warranty

This product is warranted for 1 year from the purchase date.

1.In accordance with the instructions to use, in the warranty period, it it' s not work, factory will offer free maintenance. (Please submit the goods repair and warranty card)

2.Even in the warranty period, if it's damaged caused by the following circumstances..

2.1 Malfunctions and damage caused by incorrect use

2.2Due to their own improper repair or transformation of the failure and damage

2.3After the purchase due to transport, drop, etc. arising from the failure and damage

2.4Occurrence of fire, earthquake, flood, thunder and other weather variations, abnormal voltage generated by the failure and damage

3.The warranty card will not be valid once it' s separately from the user manual, please keep it safely.

Product Warranty Card

Order Number			
Purchase Name		Phone Number	
Purchase Date		Email	
Purchase Address			
Maintenance Record			
Order Number	Maintenance Part	Maintenance Result	Maintenance Person

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